

Panel PC with Touch Function and Pre-Installed Software for Control of the EIZO Large Monitor Manager

Easy to Operate

Due to the high resolution and improved touch reaction of a modern Panel PC the CID1300P provides a fast and intuitive selection of video sources directly by the Large Monitor Manager view management. In addition, a selection via the LMM Control software is possible, which is available as an option and is installed on the Large Monitor Manager.



CID1300P

Simple Integration in the Video Management System

The introduction of simple input screens makes specialized IT administrator knowledge unnecessary. The installation effort is reduced to two steps: integration into the local network and entry of the video management system's IP address.


Additional Security Measures

Providing security measures such as service passwords, blocked network ports, and restricted system access, the CID1300P is designed for robust and secure use in the OR environment.


Specifications

Housing	Material	Aluminum powder-coated, anti-bacterial
	Color	White
Panel	Size	13.3"
	Native Resolution	1920 x 1080
	Viewable Image Size (H x V)	293 x 165 mm
	Viewing Angles	178°/178°
	Brightness	360 cd/m ²
USB	Downstream	USB 5Gb/s: Type A, USB 2.0: Type A (for connection of keyboard, mouse or updates)
Interfaces	Video Output	HDMI Type A
	Audio	Audio Jack (line-out/mic in)
	Communication	Ethernet (RJ-45): 10/100/1000 Mbit/s (Connection to video management network), RS232/422/485 galvanic isolation up to 4 kV
HW / SW Components	Operating System	Windows 10 IoT Enterprise 64-bit
	IT Components	Intel Atom x5-E3930 processor, SO-DIMM 4GB DDR3L RAM 1866 MHz, SSD via internal M.2 2280 slot (128 GB)
	Touchscreen	Multi-touch (projective capacitive), glove compatible, protective glass (7H) and anti-glare coating
Power	Power Input	AC 100-240V: 47-63Hz medical AC adapter (PSE marking)
	Power Consumption	65 W
Physical Specifications	Net Weight	3.2 kg
	Hole Spacing (VESA Standard)	VESA 100 x 100, M4 x 6 mm
Environmental Requirements		Front: IP65 Entire device: IP21
Certificates and Standards (Please contact EIZO for the latest information)		CE, EN55032, EN55035, EN62368-1, EN IEC63000, IEC 60601-1:2012 & EN 60601-1:2013, IEC 60601-1-2:2014 & EN 60601-1-2:2015, ANSI/AAMI ES60601-1:2005(R)2012, CSA C22.2 No. 60601-1:2014, FCC, REACH, RoHS
Supplied Accessories (May vary by country. Please contact EIZO for details)		AC adapter, screw set, Utility Disc (PDF Instructions for Use)

Dimensions (mm)

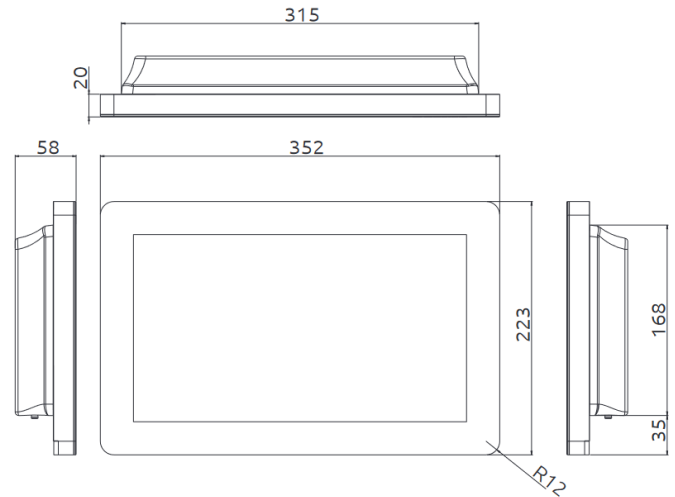

Connectors

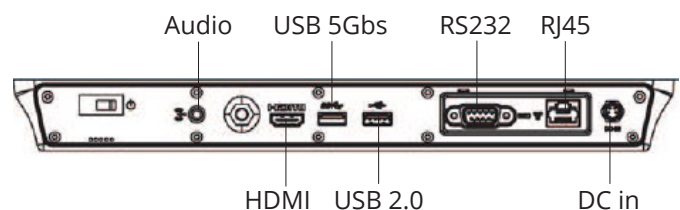

You will find your EIZO contact partner at:
<https://www.eizo-or.com/contact>

EIZO, the EIZO logo, and Curator are registered trademarks of EIZO Corporation in Japan and other countries. The terms HDMI and HDMI High-Definition Multimedia Interface, and the HDMI logo are trademarks or registered trademarks of HDMI Licensing, LLC in the United States and other countries. All other company names, product names, and logos are trademarks or registered trademarks of their respective owners. Specifications are subject to change without notice.